

O2k-Manual: Oxia



Mitochondrial Physiology Network 26.14(02):1-8 (2022)

Version 02: 2022-10-31 ©2022 Oroboros

Updates: https://wiki.orooboros.at/index.php/MiPNet26.14_Oxia

Oxia – HyperOxia to Hypoxia

Sabine Schmitt¹, Alexander Merth²,
Michael Walter-Vracevic², Erich Gnaiger¹

¹Oroboros Instruments GmbH

High-Resolution Respirometry
Schoepfstrasse 18, 6020 Innsbruck, Austria
Email: instruments@orooboros.at
www.orooboros.at

²WGT-Elektronik GmbH & Co KG

Rettenbergstraße 30a, A-6114 Kolsass, Austria




Contents

1.	Intended use.....	1
2.	Safety.....	2
3.	Components.....	3
4.	Specifications.....	4
4.1.	Technical specifications.....	4
4.2.	Ambient conditions	5
4.3.	Safety specifications.....	5
5.	Operating instructions	5
5.1.	Assembly	5
5.2.	Production and withdrawal of O ₂ and H ₂ gas.....	6
5.3.	Setting O ₂ concentrations in the O2k-chamber using O ₂ or H ₂ gas.....	7
5.4.	Storage.....	8
6.	Troubleshooting	8
7.	Author contributions and acknowledgements	8

1. Intended use

The Oxia generates gaseous oxygen and hydrogen by electrolysis of water using a proton exchange membrane PEM. O₂ and H₂ gas are used to control the O₂ regime in the Oroboros O2k. Low oxygen concentrations (<50 μM) mimic tissue normoxia or hypoxia. Hyperoxic conditions above air saturation (250-600 μM O₂) are routinely used for high-resolution respirometry of permeabilized muscle fibers or to induce oxidative stress in cells and mitochondrial preparations.

2. Safety

- Before operating the Oxia, read the user manual. Inappropriate handling can cause the formation of combustive oxyhydrogen.
- Flammable gases, Category 1. Keep the Oxia away from heat, hot surfaces, sparks, open flames, and other ignition sources. 
- The produced gas is for laboratory use only.
- Do not leave operating unattended.
- Oxia must not be used outdoor.
- Ensure easy access to power socket enabling the external power supply to be quickly and safely unplugged in case of malfunction.
- Switch off after operation and for overnight storage.
- Do not use grease on the valves.
- Do not operate the Oxia in a space smaller than 9 m³, such as a small cabinet.
- Unless there is a continuous air exchange, exchange air in the room for 5 min after 20 cycles of gas production.
- Do not cover the Oxia during operation.
- The Oxia must only be operated with the external power supply supplied or recommended by Oroboros Instruments.
- The Oxia-Base must not be exposed to liquids (IP Code 21).
- Use separate Gas syringes for withdrawal of O₂ or H₂ gas, respectively.
- Do not switch on the Oxia without the Oxia-Cell being connected to the Oxia-Base.
- Do not operate the Oxia without water in the H₂O-chamber.
- Do not use any liquid besides deionized or distilled H₂O to operate the Oxia.
- Do not disassemble the Oxia-Cell by loosening the hex head screws.
- Any servicing or repair must be done by the manufacturer. In case of any problems please contact Oroboros Instruments (<https://www.orooboros.at/index.php/o2k-technical-support/>).



Flammable gases, Category 1.




General caution




Direct current

3. Components

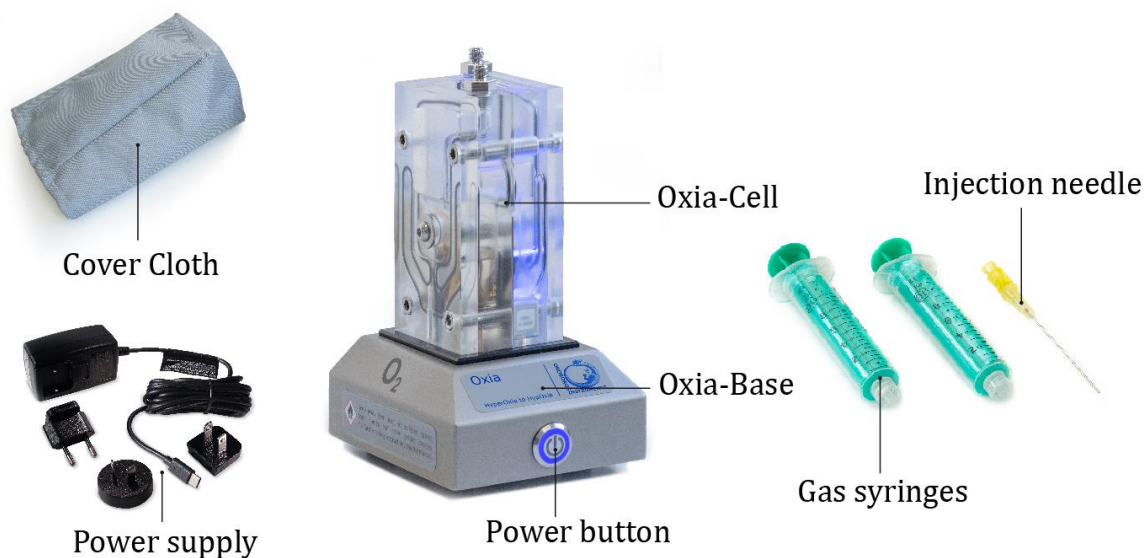


Oxia-Base

(1) Power button on the front; (2) electrical contact pins; (3) LEDs for gas chamber illumination; (4) stainless steel assembly guide; red circles: sensors to detect filling of the O₂- and H₂-chambers on top; not shown: USB-C socket (on the rear).

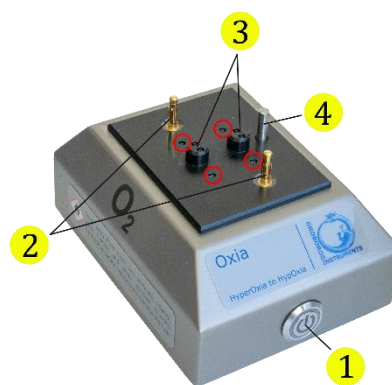

Oxia-Cell

Proton exchange membrane (PEM) electrolysis cell and (1) separate H₂O-chambers A (left) and B (right), which are connected to the (2) O₂ or (3) H₂-chamber, respectively. (4) The Luer Lock connections with spring valves allow for gas withdrawal from the O₂- and H₂-chamber by the Gas syringes.

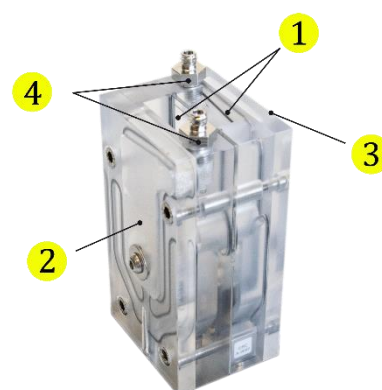

Do not adjust the Luer Lock connectors on the Oxia-Cell with a wrench or other tool as this may crack the cell.

Cloth Cover	to protect the Oxia-Cell from contamination during storage.
Gas syringes	10 mL syringes with Luer Lock connection for gas withdrawal from the O ₂ - and H ₂ -chamber and subsequent gas injection into the O ₂ k-chambers.
O₂k-Injection needles:	with spacers to obtain the correct insertion length for gas injection with the Gas syringes into the O ₂ k-chambers.
Voltage supply	The Oxia is plugged to a power socket via an external power supply.

Use only the power supply provided by the manufacturer.

Do not use a USB-A to USB-C adapter.

4. Specifications

4.1. Technical specifications

Voltage supply:	USB-C Power supply with 1.5 A, 5 V, max. 7.5 W
Dimensions:	138x108x206 mm
Weight:	1.5 kg
Gas production rate:	O ₂ : 7 mL/min H ₂ : 14 mL/min
Volume of H₂O-chambers:	37 mL
Volume of O₂-chamber:	32 mL
Volume of H₂-chamber:	32 mL
External Power Supply:	Voltage: AC 100 – 240 V Frequency: 50 – 60 Hz Current: 0.45 or 0.6 A Overvoltage category II Pollution degree 2

4.2. Ambient conditions

Operating temperature:	5-40 °C
Maximum height above sea level:	2000 m
Maximum relative air humidity:	80 % at 31 °C, 50 % at 40 °C; linear temperature dependence
Tolerance of voltage supply:	max ±10 %

4.3. Safety specifications

- Optical sensors for the O₂- and H₂-chamber stop electrolysis and thereby gas production as soon as the water is completely replaced by gas, indicated by illumination (blue) of the respective gas chamber.
- If gas production is not automatically stopped upon filling of one of the gas chambers, the excessive gas escapes to the H₂O-chamber. The H₂O-chambers are open to allow the gas to escape. The separator between the two H₂O-chambers prevents mixing of H₂ and O₂. Thus, there is no risk of oxyhydrogen combustion.
- The minimum cubature of the room (>9 m³) and frequent air exchange prevent the formation of oxyhydrogen, as long as the volume fraction of hydrogen is kept below 4 % at atmospheric pressure.

5. Operating instructions

5.1. Assembly



Ensure the Oxia power is switched off (power button is not illuminated) when the Oxia-Cell is not connected to the Oxia-Base.


1. Both H₂O-chambers and both gas chambers of the Oxia-Cell must be empty before assembly.
 - a. To remove residual water from the H₂O-, O₂-, and H₂-chambers, switch off the Oxia and unplug the Oxia-Cell from the Oxia-Base.
 - b. Remove water from the H₂O-chambers by turning the Oxia-Cell upside down.
 - c. Withdraw residual water from the O₂- and H₂-chambers with a gas syringe via the Luer Lock connection with the Oxia-Cell upside down.
2. Plug the Oxia-Cell onto the Oxia-Base. The stainless-steel assembly guide on the Oxia-Base must line up with the corresponding hole on the bottom of the Oxia-Cell to ensure correct alignment. The Oxia-Cell is flush against the Oxia-Base if inserted correctly.

5.2. Production and withdrawal of O₂ and H₂ gas



Before switching on the Oxia, ensure that the Oxia-Cell is correctly connected to the Oxia-Base.

1. Fill both H₂O-chambers with deionized or distilled H₂O up to the fill mark (1).
2. Screw an empty gas syringe with the piston fully inserted onto the Luer connector of the O₂-chamber and withdraw air from the O₂-chamber by suction (2). Thereby, the water is sucked from the corresponding H₂O-chamber into the O₂-chamber. Repeat until the water level is right beneath the spring valve, as seen from the side (3) or on top (4).
3. Repeat step 2 for the H₂-chamber.

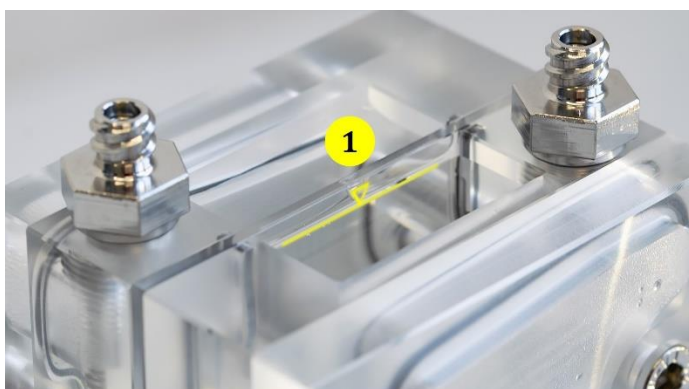

4. Connect the Oxia (USB-C socket on the rear of the Oxia-Base) to a power socket via the external power supply.
5. Start gas generation by pressing the power button on the front of the Oxia-Base. The power button is illuminated in blue indicating that the Oxia is on.
6. Let the Oxia on until the H₂-chamber is filled with gas or until enough O₂ or H₂ gas is produced for your experiment. Water electrolysis and thereby gas production stops automatically when the O₂- or H₂-chamber is filled with gas. The blue illumination is automatically switched on in the respective chamber.
7. If only one type of gas (O₂ or H₂) is needed, the other gas can be continuously released by screwing an open syringe onto the respective Luer Lock connector.
8. It is recommended to withdraw the gas right before injection to the O₂k-chamber as it mixes with air when stored in the open gas syringe. The gas can be stored for up to one hour in the gas chamber of the Oxia-Cell.


- For gas withdrawal, screw a Gas syringe to the Luer Lock connector of the Oxia-Cell. Pull out the amount of gas needed for the experiment. Screw the O2k-Injection needle immediately onto the Gas syringe. Use only the provided needle and spacer to obtain the correct insertion length for gas injection into the O2k-chambers.
- One filling of the H₂O-chambers is sufficient for about 25 cycles.
- Before refilling a H₂O-chamber the respective gas chamber must be filled completely with gas. Otherwise H₂O may spill over upon gas production.

5.3. Setting O₂ concentrations in the O2k-chamber using O₂ or H₂ gas

- These instructions apply to the 2 mL O2k-chamber.
- Lift the stopper of the O2k-chamber into the open position using the Stopper-Spacer. Gas injection into aqueous phase must be strictly avoided.
- Insert the injection needle into the titration port and gently inject the gas.
- Remove the injection needle.
- When the targeted O₂ concentration is approached, gently insert the stopper fully to close the O2k-chamber.

Further details:

https://wiki.oroboros.at/index.php/Setting_the_oxygen_concentration

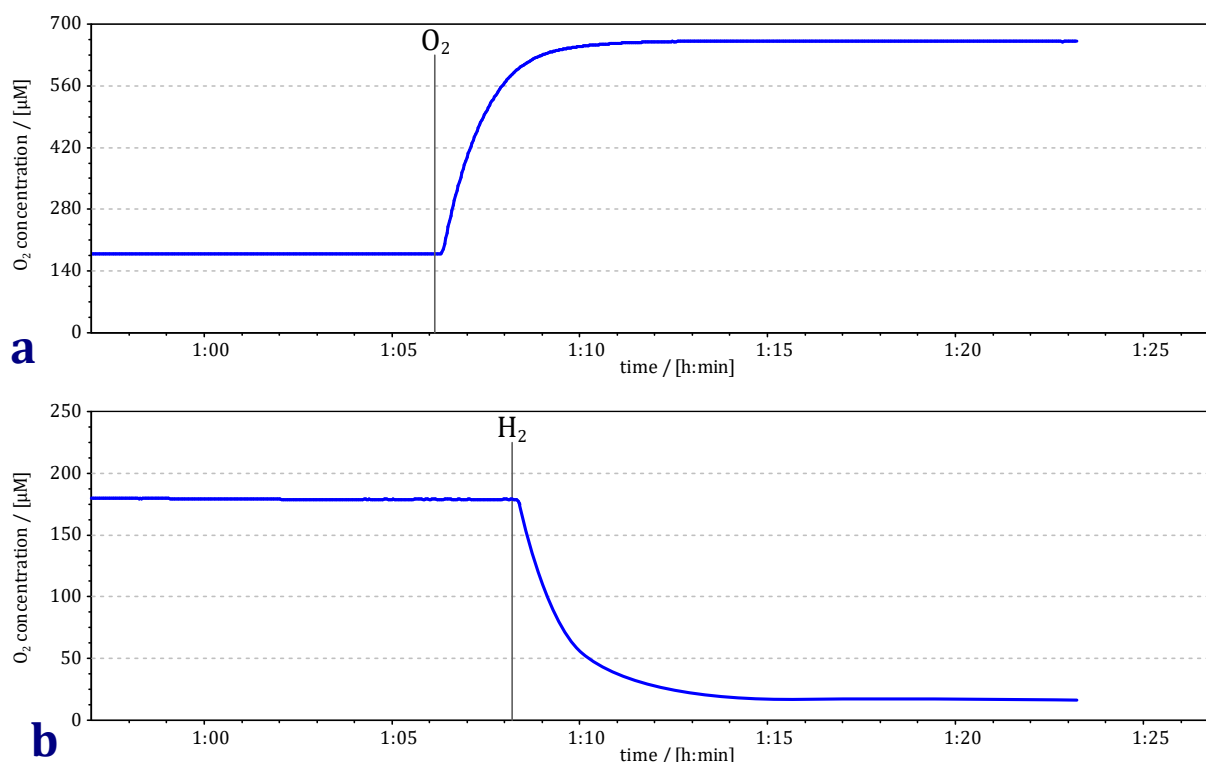

Traces of O₂ concentration in the O2k-chamber upon injection of (a) O₂ and (b) H₂ gas: The O2k-chambers were not closed to illustrate the time courses of the O₂ regime after injections of (a) 2 mL O₂ or (b) 8 mL H₂. 2 mL experimental O2k-chamber volume with MiR05 at 37 °C.

5.4. Storage

- Remove residual water from the H₂O-, O₂-, and H₂-chambers as described in 5.1, Step 1.
- Cover the Oxia-Cell with the Cloth Cover to protected from contamination. There is no need to clean the Oxia-Cell internally. Never use an alcohol-based cleaner or any strong acidic or alkaline solutions as this can damage the material of the Oxia-Cell.

6. Troubleshooting

- Air is sucked into the O₂- or H₂-chamber if the H₂O-chambers are empty during gas withdrawal. ⇒ Switch off the Oxia. Unplug the Oxia-Cell and repeat steps of section 5.1 and steps 1 and 2 of section 5.2.
- Water spills over during gas generation if too much water is added to the H₂O-chambers. ⇒ Switch off the Oxia. Unplug the Oxia-Cell and repeat steps of section 5.1 and steps 1 and 2 of section 5.2.
- Water is sucked into the Luer Lock connector. ⇒ Let the Oxia run for about 5 min. Carefully suction gas thereby removing residual water from the Luer Lock connector. Repeat if necessary.
- Power button flashes. ⇒ Ensure that the Oxia-Cell is correctly plugged onto the Oxia-Base (see section 5.1, step 2) and that the H₂O-chamber contains enough water (repeat steps from section 5.1).
- Power button is still flashing although the Oxia-Cell is correctly inserted. ⇒ Contact Oroboros Technical Support.
- Gas escapes to the H₂O-chambers. ⇒ Contact Oroboros Technical Support
- Gas production does not stop automatically when the gas chambers are filled with gas. ⇒ Contact Oroboros Technical Support.

Technical Support is provided by:

Oroboros Instruments

High-Resolution Respirometry

Schoepfstrasse 18, 6020 Innsbruck, Austria

<https://www.orooboros.at/index.php/support/>

7. Author contributions and acknowledgements

Gnaiger E was responsible for the concept of the project. Walter-Vracevic M, and Merth A were responsible for electronic and mechanical development of the Oxia. Schmitt S and Gnaiger E prepared the manual and all coauthors contributed to the final version.

We thank Lisa Tindle-Solomon for proofreading and Paolo Cocco for graphics.